



Update _2_0_181A

Clinical Reminders

VA-COVID-19 POC ANTIGEN RESULTS NOTE UPDATE

Install Guide

March 2024

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION	1
INSTALL DETAILS	8
PRE-INSTALL	9
INSTALL EXAMPLE	10
POST INSTALLATION	14

Introduction


DESCRIPTION:

This update contains three components:

1. An update to the COVID-19 POC ANTIGEN RESULTS NOTE template.

After feedback from the field with concerns that the options in the COVID-19 POC Antigen testing templates did not have health factors, the following changes have been made:

BEFORE:

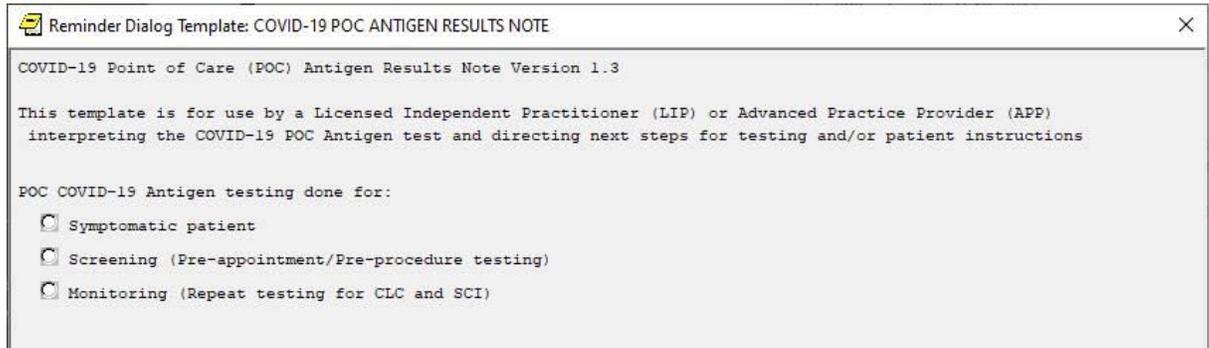

Reminder Dialog Template: COVID-19 POC ANTIGEN RESULTS NOTE


COVID-19 Point of Care (POC) Antigen Results Note Version 1.3

This template is for use by a Licensed Independent Practitioner (LIP) or Advanced Practice Provider (APP) interpreting the COVID-19 POC Antigen test and directing next steps for testing and/or patient instructions

POC COVID-19 Antigen testing done for:

- Symptomatic patient
- Screening (Pre-appointment/Pre-procedure testing)
- Monitoring (Repeat testing for CLC and SCI)

EACH OPTION EXPANDS TO POSITIVE AND NEGATIVE BUT NO HEALTH FACTORS WERE FILING FOR THE RESULTS:

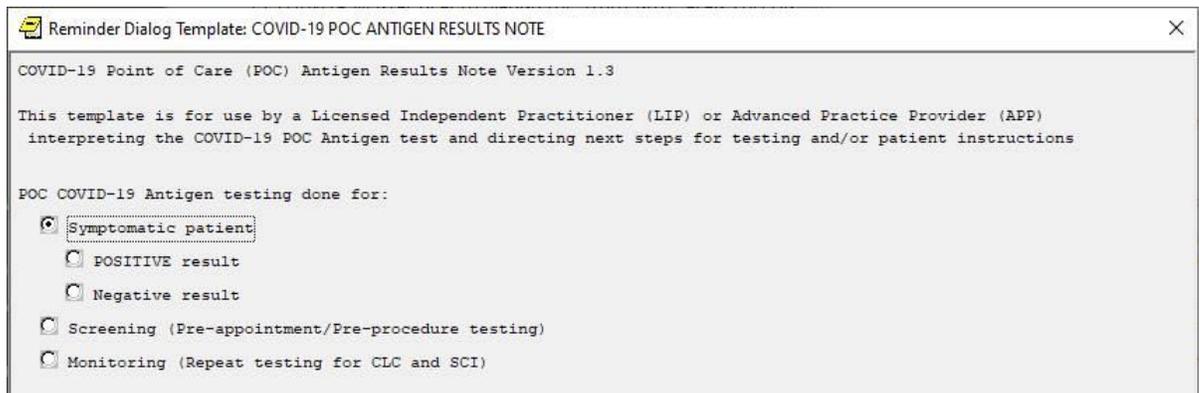

Reminder Dialog Template: COVID-19 POC ANTIGEN RESULTS NOTE

COVID-19 Point of Care (POC) Antigen Results Note Version 1.3

This template is for use by a Licensed Independent Practitioner (LIP) or Advanced Practice Provider (APP) interpreting the COVID-19 POC Antigen test and directing next steps for testing and/or patient instructions

POC COVID-19 Antigen testing done for:

- Symptomatic patient
- POSITIVE result
- Negative result
- Screening (Pre-appointment/Pre-procedure testing)
- Monitoring (Repeat testing for CLC and SCI)

The positive response advice in this template did not match the other COVID-19 templates

Before:

++Patient notified this is considered definitive testing.++

Patient with positive COVID-19 antigen result was advised to:

- Stay at home except to get medical care that cannot be provided by a telephone or video visit.
- Get rest and stay hydrated.
- Call ahead before visiting your doctor and tell them that you have or may have COVID-19.
- Wear a facemask if you are around others.
- Cover your coughs and sneezes.
- Clean your hands often.
- Keep a six-foot distance from others in your home, ideally in a separate room. Isolate in a private room with a separate bathroom, if available.
- Avoid sharing household items.
- Clean and disinfect high touch surfaces daily, ideally with a product that kills cold and flu viruses.
- If symptoms worsen, call 911.
- Contact VA nurse telephone hotline for any new or concerning symptoms at (XXX) XXX-XXXX.

Take over the counter medications for cough or fever per manufacturer recommendations.

Negative monitoring and screening before:

POC COVID-19 Antigen testing done for:

- Symptomatic patient
- Screening (Pre-appointment/Pre-procedure testing)
 - POSITIVE result
 - Negative result

++Patient notified the test was negative and no additional testing is needed++


To prevent getting sick, patient advised to:

- Continue physical distancing from others.
- Wear a mask.
- Clean hands.

Updates to the template:

- Version # updated.
- Added options to document test name, expiration and lot number with health factors for tracking purposes if needed.
- Added health factors to the positive and negative response for each of the three options.
 - o These health factors will be in the banner information but will not change the COVID-19 Banner status.
- Updated positive test advice to match the other national COVID-19 templates.
- Simplified negative advice for screening and monitoring.

After:

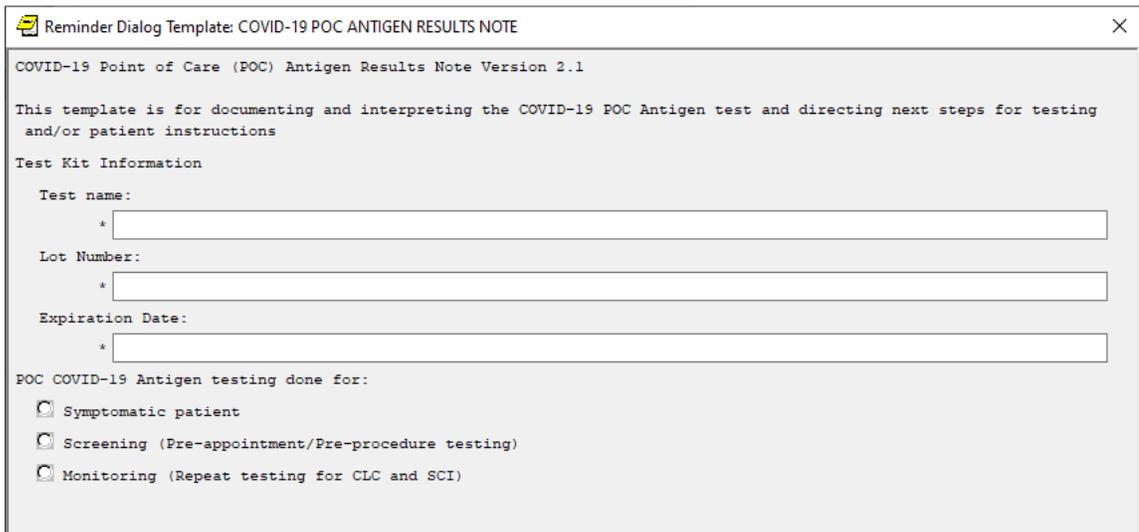

Reminder Dialog Template: COVID-19 POC ANTIGEN RESULTS NOTE

COVID-19 Point of Care (POC) Antigen Results Note Version 2.1

This template is for documenting and interpreting the COVID-19 POC Antigen test and directing next steps for testing and/or patient instructions

Test Kit Information

Test name:
*

Lot Number:
*

Expiration Date:
*


POC COVID-19 Antigen testing done for:

Symptomatic patient

Screening (Pre-appointment/Pre-procedure testing)

Monitoring (Repeat testing for CLC and SCI)

Positive results in any section and negative symptomatic display this advice which aligns with the COVID-19 discharge and other templates:

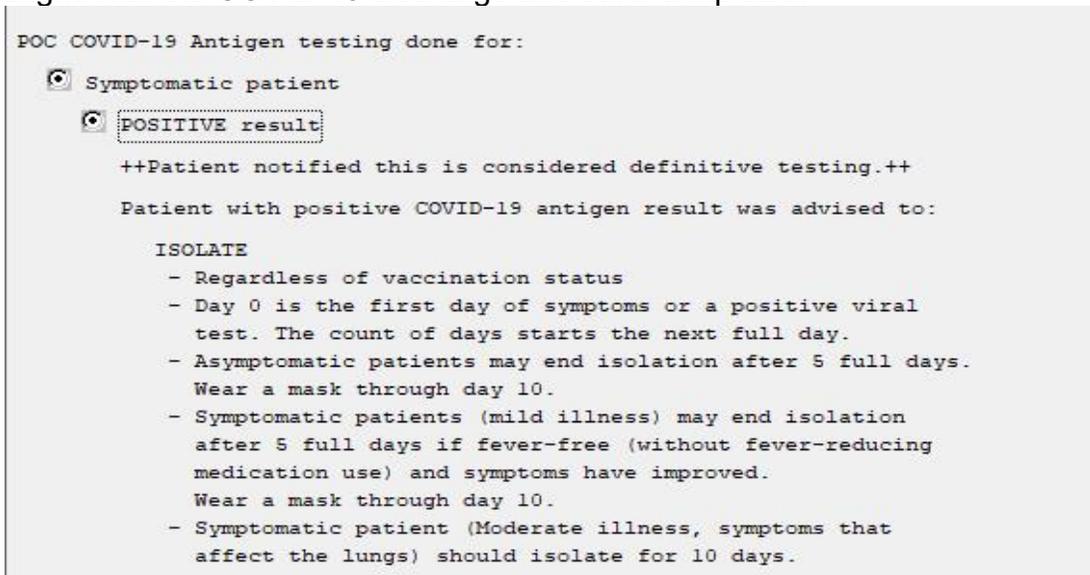

POC COVID-19 Antigen testing done for:

Symptomatic patient

POSITIVE result

++Patient notified this is considered definitive testing.++

Patient with positive COVID-19 antigen result was advised to:

ISOLATE

- Regardless of vaccination status
- Day 0 is the first day of symptoms or a positive viral test. The count of days starts the next full day.
- Asymptomatic patients may end isolation after 5 full days. Wear a mask through day 10.
- Symptomatic patients (mild illness) may end isolation after 5 full days if fever-free (without fever-reducing medication use) and symptoms have improved. Wear a mask through day 10.
- Symptomatic patient (Moderate illness, symptoms that affect the lungs) should isolate for 10 days.

Negative screening or monitoring display this advice:

POC COVID-19 Antigen testing done for:

Symptomatic patient
 Screening (Pre-appointment/Pre-procedure testing)
 POSITIVE result
 Negative result

++Patient notified the test was negative and no additional testing is needed++

To reduce likelihood of getting sick, patient advised to:

- Wear a mask around others who might be sick
- Avoid crowded places
- Clean hands

- Reminder term VA-COVID-19 CLINICAL INFO updated to include the POC result health factors listed so they will display when you click on the banner or review the health summary.

VA-COVID-19 DIAGNOSTIC ANTIGEN NEGATIVE, PCR NEEDED
 VA-COVID-19 DIAGNOSTIC ANTIGEN POSITIVE
 VA-COVID-19 SCREENING ANTIGEN NEGATIVE
 VA-COVID-19 SCREENING ANTIGEN POSITIVE, PCR NEEDED
 VA-COVID-19 MONITORING ANTIGEN NEGATIVE
 VA-COVID-19 MONITORING ANTIGEN POSITIVE, PCR NEEDED

COVID-19 Prior Positive Outside Test: 1/1/2023

The patient had a positive PCR for COVID-19 outside this VA facility on 1/1/2023.

Any non-immunocompromised patient with a prior positive PCR for COVID-19 will have a CPRS banner that displays as "Prior Positive" when 14 days has passed from the most recent positive result unless a clinician has entered that the patient is "Clinically Resolved".

LAB RESULTS:

Date	Test Name	Result
09/22/2022@06:03:23	COVID-19	Detected

No Ab lab results

OUTSIDE COVID-19 LABS:

Date	Health Factor Name	Comment	Location
01/01/2023	VA-COVID-19 PCR LAB OUTSIDE POSITIVE	SELF REPORT	
08/26/2021	VA-COVID-19 UNK OUTSIDE LAB TEST NEG	SELF REPORT	
06/02/2020	VA-COVID-19 PCR LAB OUTSIDE POSITIVE	SELF REPORT	

OTHER CLINICAL INFORMATION:

Date	Health Factor Name
02/06/2024	VA-COVID-19 MONITORING ANTIGEN NEGATIVE
05/25/2020	VA-COVID-19 RESOLVED

COVID-19 Problem List:
None

Date	Immunization	Series
09/21/2023	COVID-19 (PFIZER), MRNA, LNP-S, PF, TRIS-SUC	1
09/01/2023	COVID-19 (MODERNA), MRNA, LNP-S, PF, 50 MCG/	1
09/22/2022	COVID-19 (PFIZER), MRNA, LNP-S, BIVALENT, PF	1
08/01/2022	COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0	3
04/01/2022	COVID-19 (JANSSEN), VECTOR-NR, RS-AD26, PF,	1
04/01/2022	COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0	2

Print Close

3. A new health summary is included that duplicates the information found in the COVID-19 Banner display.
 - a. It will show each true finding and then a summary at the end that matched that click on the banner
 - b. The new health summary should be added to the reports tab in CPRS:

COVID-19 Prior Positive Outside Test: 1/1/2023

Available Reports

- > Clinical Reports
- > Health Summary
 - Adhoc Report
 - Covid-19 Information
 - Va Remote Teledermatology
 - Va Remote Tele-Eye Image
- > HDR Reports
- > Dept. of Defense Reports
 - Imaging (local only)
 - Graphing (local only)
 - Lab Status
 - Blood Bank Report
- > Anatomic Pathology
- Dietetics Profile
- Nutritional Assessment
- Vitals Cumulative
- Procedures (local only)
- Daily Order Summary
- Order Summary for a Date Range
- Chart Copy Summary
- Outpatient RX Profile
- Med Admin Log (BCMA)
- Med Admin History (BCMA)
- Surgery (local only)
- > Event Capture

- Dietetics Profile
- Nutritional Assessment
- Vitals Cumulative
- Procedures (local only)
- Daily Order Summary
- Order Summary for a Date Range
- Chart Copy Summary
- Outpatient RX Profile
- Med Admin Log (BCMA)
- Med Admin History (BCMA)
- Surgery (local only)
- > Event Capture

- Adhoc Report
- Covid-19 Information
- Va Remote Teledermatology
- Va Remote Tele-Eye Image

Health Summary Covid-19 Information

----- CF - Reminders Findings -----

----- CPRS COVID INFORMATION -----

Information:
 Computed Finding: VA-Patient Age
 02/13/2024@13:58:51 value - 76

Reminder Term: VA-COVID-19 RECOVERED
 Health Factor: Covid-19 Resolved
 06/25/2020

Reminder Term: VA-COVID-19 PCR/AG LAB RESULTS
 Laboratory test: COVID-19; specimen: NASOPHARYNX
 09/22/2022@06:03:23 value - Detected

Reminder Term: VA-COVID-19 OUTSIDE RESULTS
 Health Factor: Covid-19 PCR Lab Outside Positive
 01/01/2023
 Comments: SELF REPORT
 06/02/2020
 Comments: SELF REPORT

Health Factor: Covid-19 Unk Outside Lab Test Neg
 08/26/2021
 Comments: SELF REPORT

Reminder Term: VA-COVID-19 CLINICAL INFO
 Health Factor: POC COVID-19 Antigen Monitoring Test Negative Result
 02/06/2024@13:57:07

Health Factor: Covid-19 Resolved
 05/25/2020

Reminder Term: VA-SARS-COV-2 VACCINES
 Immunization: COVID-19 (PFIZER), MRNA, LNP-S, PF, TRIS-SUCROSE, 30 MCG/0.3 ML
 (AGES 12+ YEARS)
 09/21/2023@09:01 series - SERIES 1

Immunization: COVID-19 (MODERNA), MRNA, LNP-S, PF, 50 MCG/0.5 ML (AGES 12+ YEARS)
 09/01/2023 series - SERIES 1

Immunization: COVID-19 (PFIZER), MRNA, LNP-S, BIVALENT, PF, 30 MCG/0.3 ML DOSE
 09/22/2022@05:24:23 series - SERIES 1
 Comments: PFR; LOT; mm/dd/yyyy

Immunization: COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0.3 ML DOSE
 08/01/2022 series - SERIES 3
 04/01/2022 series - SERIES 2
 06/30/2021 series - SERIES 2

Immunization: COVID-19 (JANSSEN), VECTOR-NR, RS-AD26, PF, 0.5 ML
 04/01/2022 series - SERIES 1
 09/05/2021 series - SERIES 1

Immunization: COVID-19 (MODERNA), MRNA, LNP-S, PF, 100 MCG/0.5ML DOSE OR 50 MCG/0.25ML DOSE
 01/02/2022 series - SERIES 1

SUMMARY OF COVID-19 INFORMATION LOCATED ABOVE

LAB RESULTS:

Date	Test Name	Result
09/22/2022@06:03:23	COVID-19	Detected
Date	Test Name	Result
No Ab lab results		

OUTSIDE COVID-19 LABS:

Date	Health Factor Name	Comment	Location
01/01/2023	VA-COVID-19 PCR LAB OUTSIDE POSITIVE	SELF REPORT	
09/26/2021	VA-COVID-19 UNK OUTSIDE LAB TEST NEG	SELF REPORT	
06/02/2020	VA-COVID-19 PCR LAB OUTSIDE POSITIVE	SELF REPORT	

OTHER CLINICAL INFORMATION:

Date	Health Factor Name
02/06/2024	VA-COVID-19 MONITORING ANTIGEN NEGATIVE
05/25/2020	VA-COVID-19 RESOLVED

COVID-19 Problem List:
None

Date	Immunization	Series
09/21/2023	COVID-19 (PFIZER), MRNA, LNP-S, PF, TRIS-SUC	1
09/01/2023	COVID-19 (MODERNA), MRNA, LNP-S, PF, 50 MCG/	1
09/22/2022	COVID-19 (PFIZER), MRNA, LNP-S, BIVALENT, PF	1
09/01/2022	COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0	3
04/01/2022	COVID-19 (JANSSEN), VECTOR-NR, RS-AD26, PF,	1
04/01/2022	COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0	2
01/02/2022	COVID-19 (MODERNA), MRNA, LNP-S, PF, 100 MCG	1
09/05/2021	COVID-19 (JANSSEN), VECTOR-NR, RS-AD26, PF,	1
06/30/2021	COVID-19 (PFIZER), MRNA, LNP-S, PF, 30 MCG/0	2

*** END ***** CONFIDENTIAL COVID-19 INFORMATION SUMMARY pg. 1 *****

Cover Sheet Problems Meds Orders Notes Consults Surgery D/C Summ Labs Reports

UPDATE_2_0_181 contains 1 Reminder Exchange entry:
UPDATE_2_0_181A VA-COVID-19 POC ANTIGEN RESULTS NOTE UPDATE

The exchange file contains the following components:

HEALTH FACTORS

VA-REMINDER UPDATES [C]
VA-UPDATE_2_0_181A
VA-COVID-19 [C]
VA-COVID-19 UNK OUTSIDE LAB TEST POS
VA-COVID-19 NON PCR AG LAB OUTSIDE POS
VA-COVID-19 MONITORING ANTIGEN NEGATIVE
VA-COVID-19 SCREENING ANTIGEN NEGATIVE
VA-COVID-19 MONITORING ANTIGEN POSITIVE, PCR NEEDED
VA-COVID-19 SCREENING ANTIGEN POSITIVE, PCR NEEDED
VA-COVID-19 DIAGNOSTIC ANTIGEN NEGATIVE, PCR NEEDED
VA-COVID-19 DIAGNOSTIC ANTIGEN POSITIVE
VA-COVID-19 NO LONGER SUSPECTED
VA-COVID-19 INCORRECT PRIOR POSITIVE PCR
VA-COVID-19 RESOLVED
VA-COVID-19 SUSPECTED
VA-COVID-19 UNK OUTSIDE LAB TEST NEG
VA-COVID-19 NON PCR AG LAB OUTSIDE NEG
VA-COVID-19 OUTSIDE PCR SPEC COLLECTED
VA-COVID-19 PCR LAB OUTSIDE NEGATIVE
VA-COVID-19 PCR LAB OUTSIDE POSITIVE
VA-COVID-19 POC ANTIGEN TEST EXPIRATION DATE
VA-COVID-19 POC ANTIGEN TEST LOT NUMBER
VA-COVID-19 POC ANTIGEN TEST NAME

REMINDER SPONSOR

NATIONAL CLINICAL REMINDER COMMITTEE
NATIONAL CENTER FOR HEALTH PROMOTION AND DISEASE PREVENTION
(NCP)
VA NATIONAL CENTER FOR HEALTH PROMOTION AND DISEASE
PREVENTION (NCP)
POPULATION HEALTH

REMINDER TAXONOMY

VA-IMMUNOCOMPROMISED FOR COVID-19 ADDL DOSE
VA-HIV INFECTION
LONG TERM SYSTEMIC STEROID
VA-PNEUMOC DZ RISK - CHEMOTHERAPY
VA-COVID-19 SNOMED CODES

REMINDER TERM

VA-REMINDER UPDATE_2_0_181A
VA-CD4 COUNT
VA-STEROID PREDN EQUIV 5-6.6MG
VA-STEROID PREDN EQUIV 6.6-9MG
VA-STEROID PREDN EQUIV 10-19MG
VA-STEROID PREDN EQUIV >=20MG
VA-CHEMO/IMMUNOSUPPRESSIVE DRUGS
VA-IMMUNOCOMPROMISED FOR COVID-19 BANNER
VA-COVID-19 PCR LAB RESULTS PRESUMP POSITIVE
VA-COVID-19 AG LAB RESULTS POSITIVE HF
VA-COVID-19 RELATED ORDERABLE ITEMS FOR BANNER DISPLAY
VA-SARS-COV-2 VACCINES
VA-COVID-19 PROB LIST DX
VA-COVID-19 CLINICAL INFO
VA-COVID-19 OUTSIDE RESULTS
VA-COVID-19 ANTIBODY LAB TEST
VA-COVID-19 PCR/AG LAB RESULTS
VA-COVID-19 NO LONGER PRESUMED
VA-COVID-19 INCORRECT PRIOR POSITIVE PCR RESULT
VA-COVID-19 OUTSIDE PCR SPECIMEN COLLECTED
VA-COVID-19 RECOVERED
VA-COVID-19 PCR LAB ORDERS
VA-COVID-19 PRESUMED

REMINDER DEFINITION

VAL-IMMUNOCOMPROMISING DRUGS/LABS/DX FOR COVID-19
VA-COVID-19 CPRS STATUS OBJECT

HEALTH SUMMARY TYPE

COVID-19 INFORMATION

HEALTH SUMMARY OBJECTS

COVID-19 BANNER INFO

TIU DOCUMENT DEFINITION

COVID-19 BANNER DISPLAY

REMINDER DIALOG

VA-COVID-19 BANNER DISPLAY OBJECT
VA-COVID-19 POC ANTIGEN RESULTS NOTE

Install Details

This update is being distributed as a web host file. The address for the host file is:
https://REDACTED/UPDATE_2_0_181A.PRD

The file will be installed using Reminder Exchange, programmer access is not required.

This update can be loaded with users on the system. Installation will take less than 5 minutes.

Pre-Install

The following elements may have local modifications, please review pre-installation so you can put back any local changed post install.

1. The reminder dialog element **VAL-COVID-19 ANTIGEN TESTING ORDER** should **have your local order for PCR testing as a finding item**.
2. The reminder dialog elements below may have local content embedded. Enabling and adding local content was an optional part of the original release.
VAL-COVID-19 ANTIGEN TESTING POS LOCAL INFO
VAL-COVID-19 ANTIGEN TESTING SYMPTOMATIC PATIENT NEG
3. The following reminder terms are included in this update. You should not overwrite any of them during installation. If you would like to do a pre-install inquiry and save for later reference in case you accidentally overwrite, you may do so.

VA-COVID-19 ANTIBODY LAB TEST
VA-COVID-19 PCR LAB ORDERS
VA-COVID-19 PCR LAB RESULTS PRESUMP POSITIVE
VA-COVID-19 PCR/AG LAB RESULTS
VA-COVID-19 RELATED ORDERABLE ITEMS FOR BANNER DISPLAY
VA-CD4 COUNT
VA-STEROID PREDN EQUIV 5-6.6MG
VA-STEROID PREDN EQUIV 6.6-9MG
VA-STEROID PREDN EQUIV 10-19MG
VA-STEROID PREDN EQUIV >=20MG
VA-CHEMO/IMMUNOSUPPRESSIVE DRUGS

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                 LMM  Load MailMan Message
CMM  Create MailMan Message           LR   List Reminder Definitions
DFE  Delete Exchange File Entry       LWH  Load Web Host File
IFE  Install Exchange File Entry       RI   Reminder Definition Inquiry
IH   Installation History              RP   Repack
Select Action: Next Screen// LWH  Load Web Host File
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

https://REDACTED/UPDATE_2_0_181A.PRD

You should see a message at the top of your screen that the file loaded successfully

Search and locate an entry titled **UPDATE_2_0_181A VA-COVID-19 POC ANTIGEN RESULTS UPDATE**

```
+Item  Entry                                     Source                                     Date Packed
313    UPDATE_2_0_181A VA-COVID-19
      POC ANTIGEN RESULTS UPDATE
314    UPDATE_2_0_182 VA-OSP SUICIDE
      SAFETY PLAN UPDATE
315    UPDATE_2_0_182 VA-OSP SUICIDE
      SAFETY PLAN UPDATE
316    UPDATE_2_0_183 VA-SCARCE
      RESOURCE ALLOCATION UPDATE
317    UPDATE_2_0_184 VA-COVID-19
      TEMPLATE UPDATES
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                 LMM  Load MailMan Message
CMM  Create MailMan Message           LR   List Reminder Definitions
DFE  Delete Exchange File Entry       LWH  Load Web Host File
IFE  Install Exchange File Entry       RI   Reminder Definition Inquiry
IH   Installation History              RP   Repack
Select Action: Next Screen// IFE  Install Exchange File Entry
Enter a list or range of numbers (1-905): 313
```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry **UPDATE_2_0_181A VA-COVID-19 POC ANTIGEN RESULTS UPDATE**. (in this example it is entry 313 it will vary by site).

The date of the exchange file should be 02/21/2024.

```
Package Version: 2.0P65

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-COVID-19 POC ANTIGEN RESULTS NOTE
  VA-COVID-19 BANNER DISPLAY OBJECT

REMINDER TERM
  VA-REMINDER UPDATE_2_0_181A

Keywords:

Components:
+      Enter ?? for more actions      >>>
IA  Install all Components           IS  Install Selected Component
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**
You will see several prompts, for all new entries you will choose **I to Install**
If you are prompted to overwrite **any other components, choose S to SKIP with one exception:**

You will overwrite ONLY one existing component, the reminder term VA-COVID-19 CLINICAL INFO:
REMINDER TERM entry named VA-COVID-19 CLINICAL INFO already exists but the packed component is different, what do you want to do?

- Select one of the following:
- C Create a new entry by copying to a new name
 - M Merge findings
 - O Overwrite the current entry
 - U Update
 - Q Quit the install
 - S Skip, do not install this entry

Enter response: **O// verwrite the current entry**
Are you sure you want to overwrite? N// **Y**

It is very important that you do not overwrite any other reminder terms that may have your local ITEMS mapped while completing the installation.

You will be prompted to install the FIRST reminder dialog component – You DO NOT need to install this item. It is for transporting the object only. *There is no harm if you do install it but it does not need to be installed.*

```
Packed reminder dialog: VA-COVID-19 BANNER DISPLAY OBJECT

Item  Seq.  Dialog Findings                                     Type  Exists
-----
1      VA-COVID-19 BANNER DISPLAY OBJECT                 element
      Finding: *NONE*

+ Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details  DT  Dialog Text      IS  Install Selected
DF  Dialog Findings  DU  Dialog Usage     QU  Quit
DS  Dialog Summary  IA  Install All

Select Action: Quit// Q
```

At the **Select Action** prompt, type Q to quit.

You will be prompted to install the reminder dialog component – You will choose to install all even though some components may have an X to indicate that they exist.

```
Packed reminder dialog: VA-COVID-19 POC ANTIGEN RESULTS NOTE [NATIONAL DIALOG]
+-----+-----+-----+-----+-----+-----+
Item  Seq. Dialog Findings                                     Type  Exists
-----+-----+-----+-----+-----+-----+
1      1  VA-COVID-19 POC ANTIGEN RESULTS NOTE                       dialog  X
2      2  5 VA-COVID-19 POC ANTIGEN RESULTS NOTE VERSION             element  X
      3  Finding: *NONE*
3      3  8 VA-GP COVID-19 ANTIGEN TEST INFO                          group    X
      4  Finding: *NONE*
4      4  8.5 VAL-COVID-19 ANTIGEN TEST TEST NAME                     element  X
      5  Finding: VA-COVID-19 POC ANTIGEN TEST NAME (HEALTH FACTOR)
5      5  COVID-19 POC ANTIGEN TEST NAME                             prompt   X
6      6  8.10 VAL-COVID-19 ANTIGEN TEST LOT NUMBER                  element  X
      7  Finding: VA-COVID-19 POC ANTIGEN TEST LOT NUMBER (HEALTH
      8  FACTOR)
7      7  COVID-19 POC ANTIGEN LOT NUMBER                            prompt   X
8      8  8.15 VAL-COVID-19 ANTIGEN TEST EXPIRATION                  element  X
+-----+-----+-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage     QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-COVID-19 POC ANTIGEN RESULTS NOTE**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```
Packed reminder dialog: VA-COVID-19 POC ANTIGEN RESULTS NOTE [NATIONAL DIALOG]
VA-COVID-19 POC ANTIGEN RESULTS NOTE (reminder dialog) installed from exchange f
+-----+-----+-----+-----+-----+-----+
Item  Seq. Dialog Findings                                     Type  Exists
-----+-----+-----+-----+-----+-----+
1      1  VA-COVID-19 POC ANTIGEN RESULTS NOTE                       dialog  X
2      2  5 VA-COVID-19 POC ANTIGEN RESULTS NOTE VERSION             element  X
      3  Finding: *NONE*
3      3  8 VA-GP COVID-19 ANTIGEN TEST INFO                          group    X
      4  Finding: *NONE*
4      4  8.5 VAL-COVID-19 ANTIGEN TEST TEST NAME                     element  X
      5  Finding: VA-COVID-19 POC ANTIGEN TEST NAME (HEALTH FACTOR)
5      5  COVID-19 POC ANTIGEN TEST NAME                             prompt   X
6      6  8.10 VAL-COVID-19 ANTIGEN TEST LOT NUMBER                  element  X
      7  Finding: VA-COVID-19 POC ANTIGEN TEST LOT NUMBER (HEALTH
      8  FACTOR)
7      7  COVID-19 POC ANTIGEN LOT NUMBER                            prompt   X
8      8  8.15 VAL-COVID-19 ANTIGEN TEST EXPIRATION                  element  X
+-----+-----+-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage     QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q
```

After completing this dialog install, you will type Q.

```
Package Version: 2.0P65

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-COVID-19 POC ANTIGEN RESULTS NOTE
  VA-COVID-19 BANNER DISPLAY OBJECT

REMINDER TERM
  VA-REMINDER UPDATE_2_0_181A

Keywords:

Components:
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
IA  Install all Components          IS  Install Selected Component
Select Action: Next Screen// C
```

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Install complete.

Post Installation

1. Confirm the COVID-19 POC ANTIGEN RESULT NOTE version displays as 2.1.
2. Review your pre-install inquiry and add your local quick order for COVID-19 PCR testing as a finding item to the reminder dialog element **VAL-COVID-19 ANTIGEN TESTING ORDER**.
3. Review your pre-install inquiry and add local content back if needed - **OPTIONAL**
 This template was built with the standard follow up instructions included. If a site would like to add guidance or information specific to that site, you can add text to the local reminder dialog element VAL-COVID-19 ANTIGEN TESTING POS LOCAL INFO for patients with a positive result and VAL-COVID-19 ANTIGEN TESTING SYMPTOMATIC PATIENT NEG for symptomatic patients with negative POC results. These dialog elements are at the end of the sections. This dialog elements will install disabled and you will need to remove the disable flag for them to be used.

4. Add the new Health Summary to your reports tab. (REQUIRED)

In VistA, using the Health Summary Coordinator's menu, choose the option CPRS Reports Tab 'Health Summary Types List' Menu ..., then choose option Edit 'Health Summary Types List' Parameters, at the system level, review the current sequence options and then add the Health Summary COVID-19 INFORMATION to a sequence as close to the top as you can.

Here is an example, the sequences will vary by site.

Select Health Summary Coordinator's Menu <TEST ACCOUNT> Option: 5
CPRS Reports Tab 'Health Summary Types List' Menu

- 1 Display 'Health Summary Types List' Defaults
- 2 Precedence of 'Health Summary Types List'
- 3 Method of compiling 'Health Summary Types List'
- 4 Edit 'Health Summary Types List' Parameters**

Select CPRS Reports Tab 'Health Summary Types List' Menu <TEST ACCOUNT> Option: **4 Edit 'Health Summary Types List' Parameters**

Allowable Health Summary Types may be set for the following:

- 2 User USR [choose from NEW PERSON]
- 3 Division DIV [choose from INSTITUTION]
- 4 System SYS [YOURSITE.MED.VA.GOV]

Enter selection: **4 System YOURSITE.MED.VA.GOV**

- Setting Allowable Health Summary Types for System:
YOURSITE.MED.VA.GOV -
Select Sequence: ?

Sequence Value

- 1 GMTS HS ADHOC OPTION
- 5 REMOTE CLINICAL DATA (1Y)
- 10 REMOTE DEMO/VISITS/PCE (3M)

Select Sequence: **2**


Are you adding 2 as a new Sequence? Yes// **YES**

Sequence: 2// 2 **<ENTER>**

Health Summary: **COVID-19 INFORMATION**

Select Sequence: **HIT ENTER TO LEAVE THIS OPTION**

5. Close and reopen CPRS, go to your reports tab and confirm that you see the new option displayed under the Health Summary in the order that you added it:

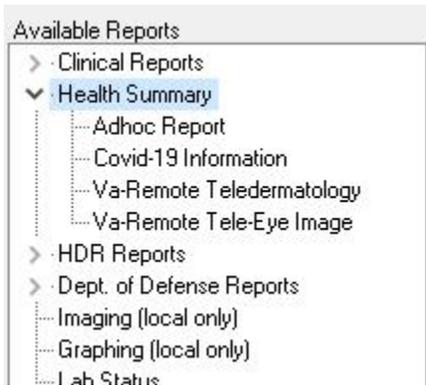

NOTE: If there are any questions or problems during the installation, please notify your local IT support to log a national help desk ticket or contact the VA Service Desk and have them submit a national ticket to:

Category: Affected Service

Affected Service: VistA – CPRS: Clinical Reminders

Assignment Group: NTL SUP Clin 2